


# P-SDTR Sensor GUI v1

User Guide

# Features


- RGB and laser projectors
- Heliotis Camera
- CMOS Camera
- Data Analysis
- Signal Generator
- Grains Finder

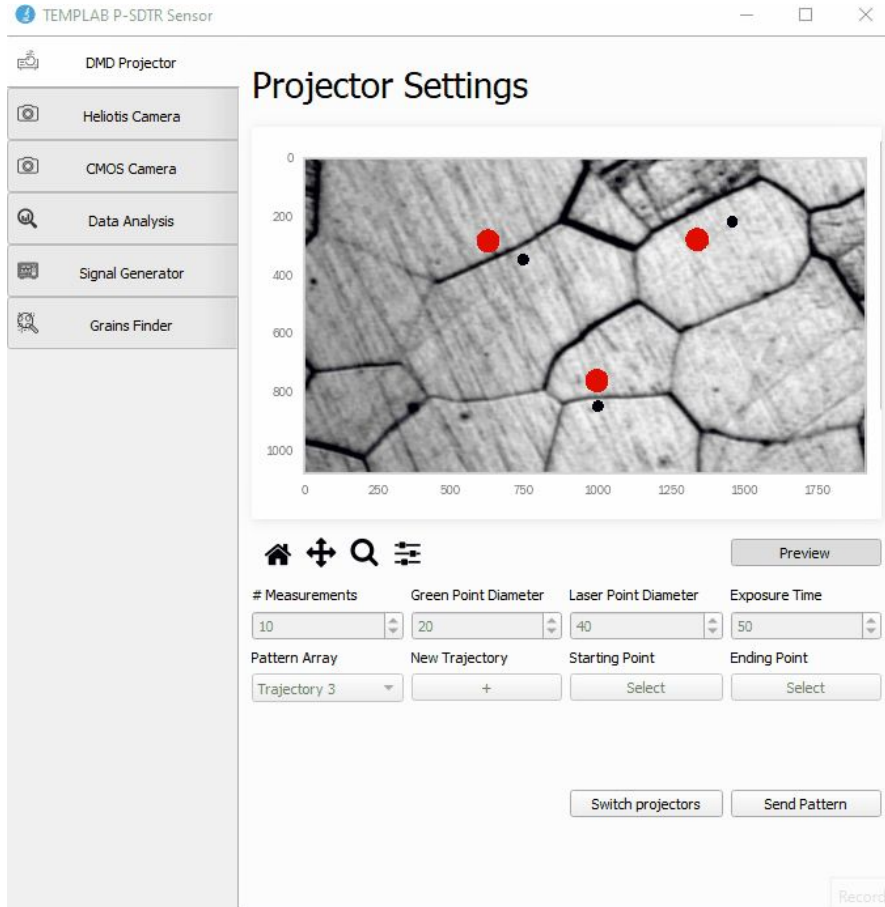

## RGB and Laser Projectors


The user configures a set of trajectories to be projected on the sample.

For each trajectory, the RGB projector will display a sequence of equally spaced points all across, and the laser projector will display a stationary point at the middle.

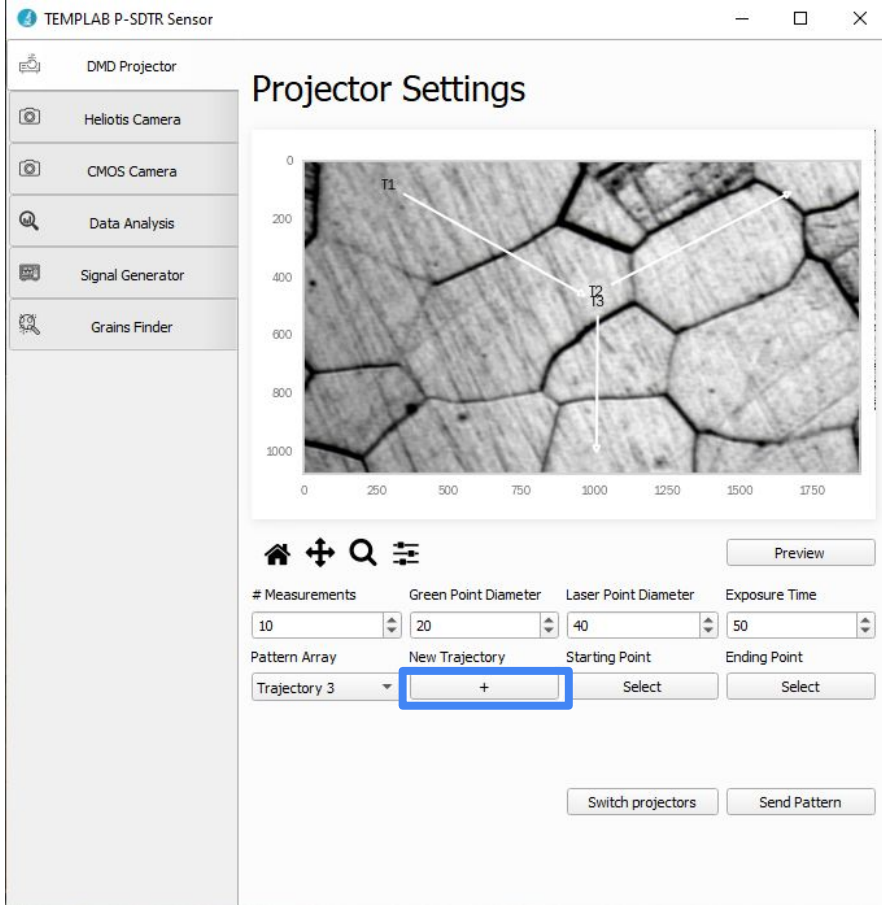

## RGB and Laser Projectors


Add a new trajectory with the “+” button.

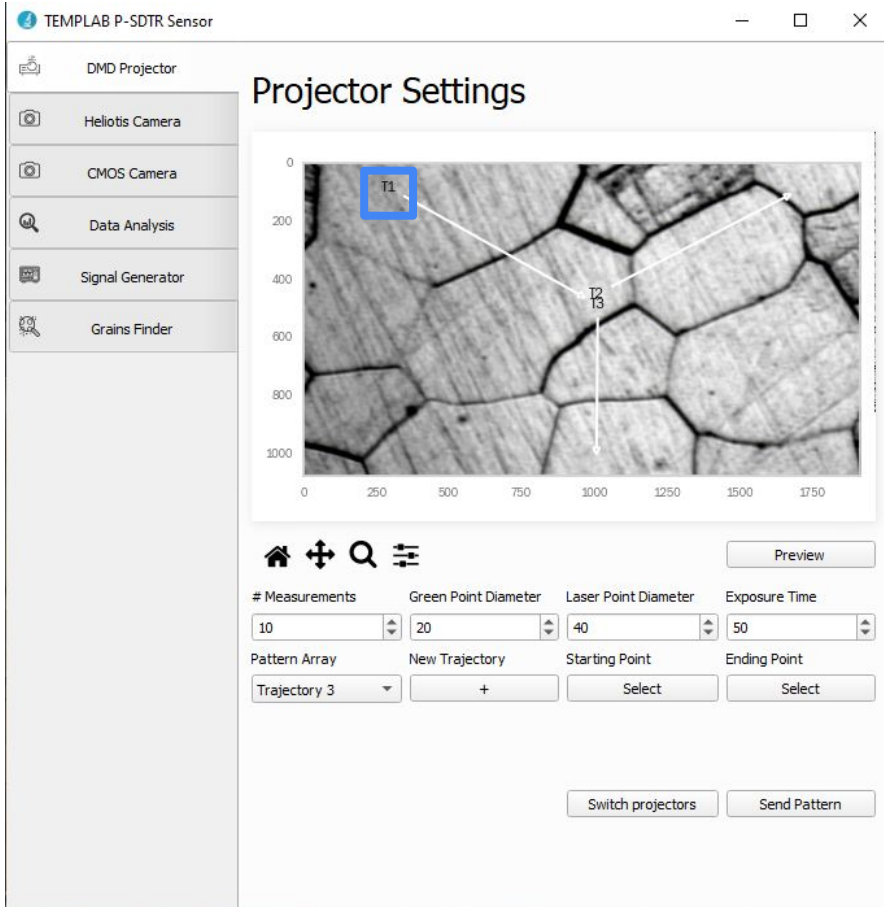

## RGB and Laser Projectors


Select the starting and ending points on the plot.

You can edit those later for any trajectory using the “Pattern Array” dropdown and the “Select” buttons.

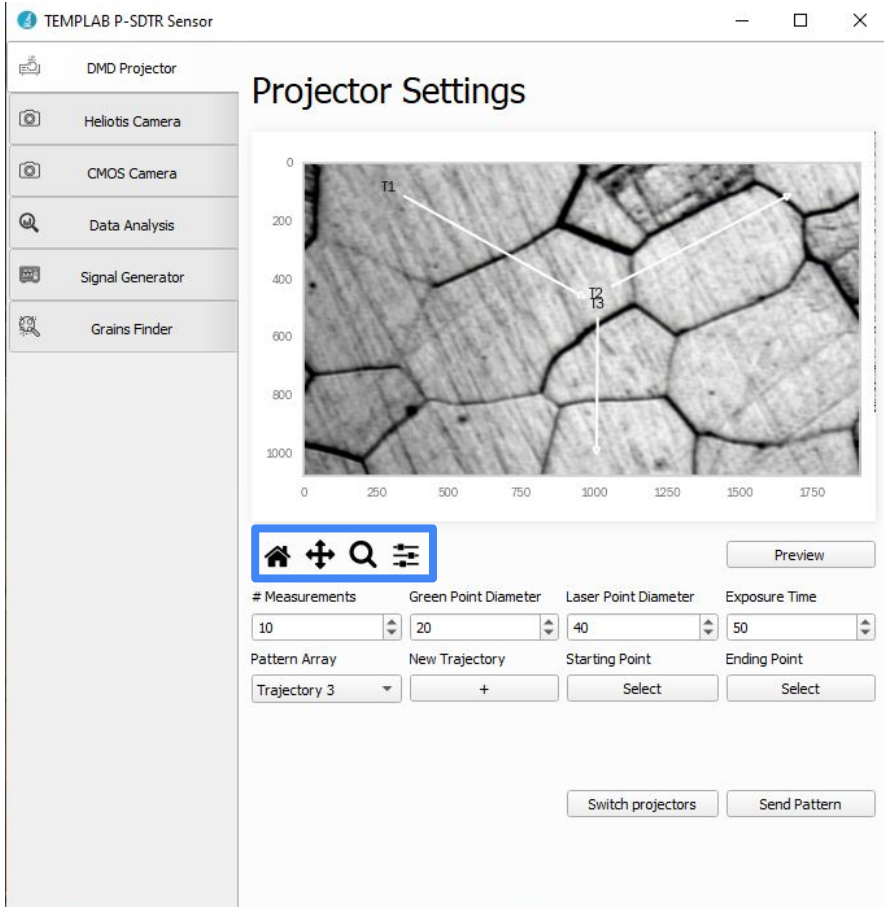


## RGB and Laser Projectors

You can interact with the plot (zoom, pan, etc) using the navigation bar.

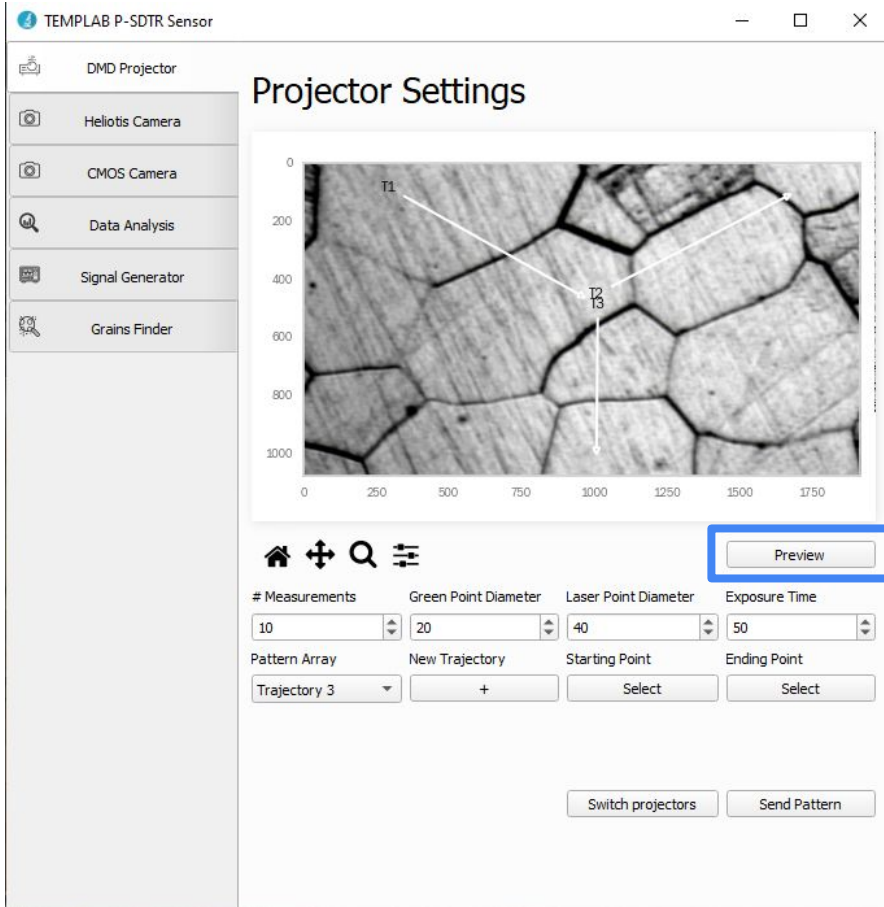

# RGB and Laser Projectors

Click “Preview” to double-check.

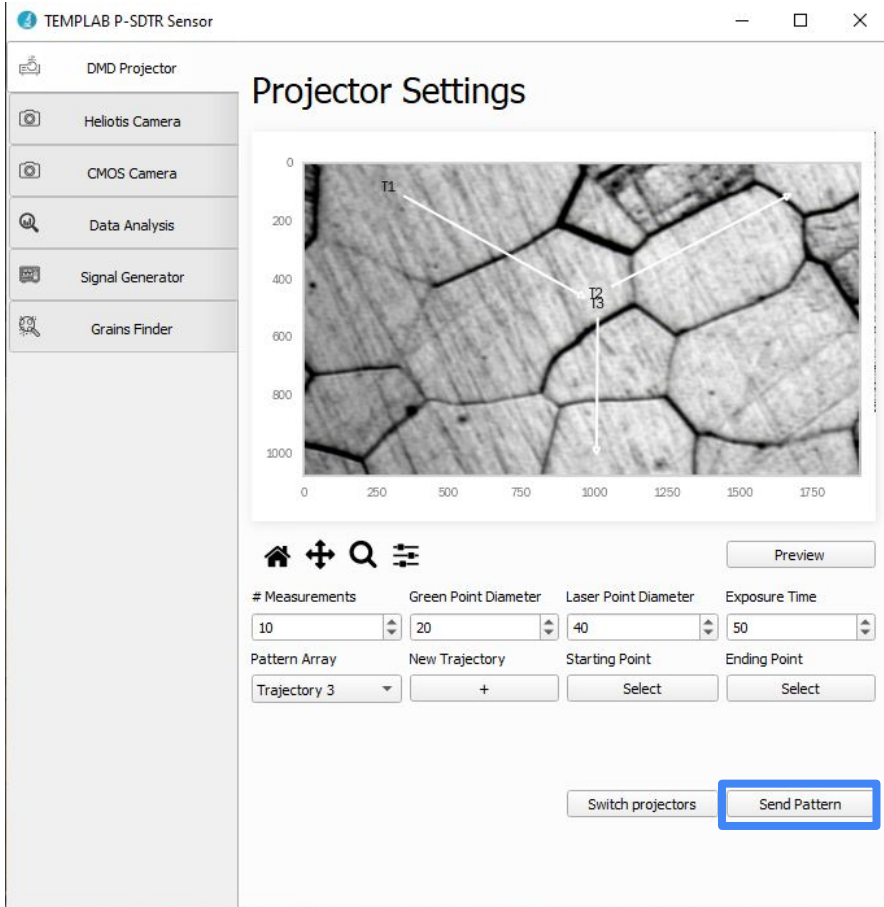


## RGB and Laser Projectors

With the devices connected, click “Send Pattern” to set up the sequence in the projectors.

If the projectors are not connected or recognized by the device, a message will appear on the Python terminal.



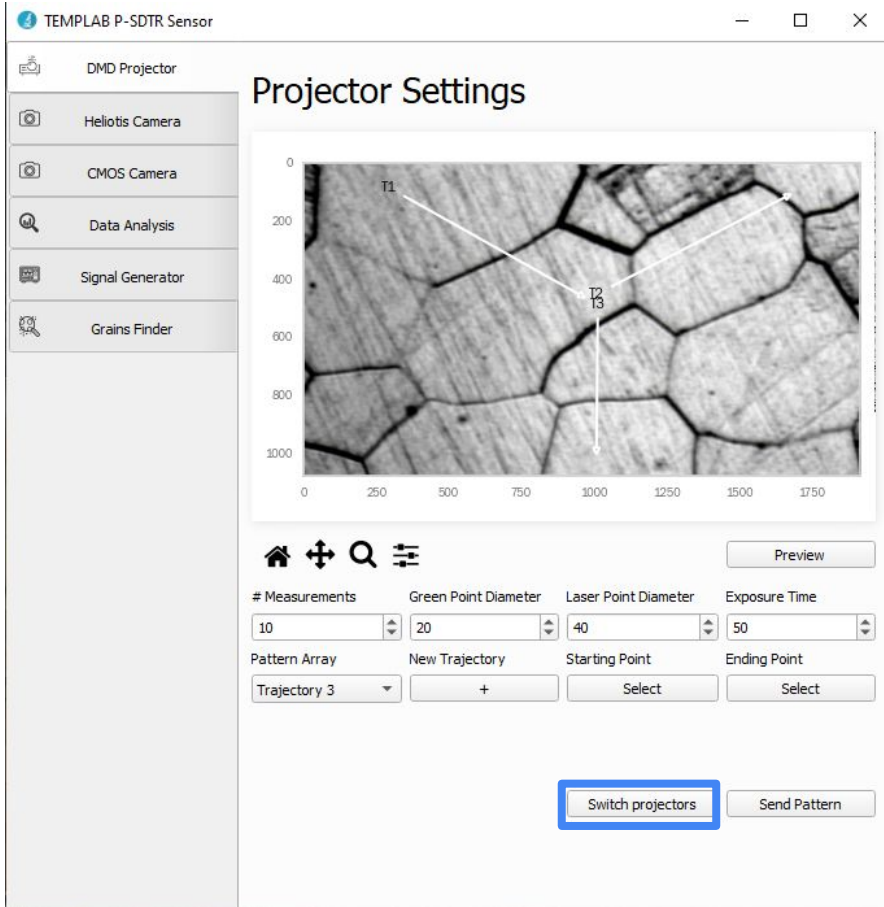

## RGB and Laser Projectors

If the order of the projectors is incorrect (laser projector showing the rgb pattern and vice versa), click on “Switch projectors” and send the pattern again.